

This Page Is Inserted by IFW Operations
and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

BioPacTM

Base Transcription

Gene

Microarray

PCR

RT-PCR

Western Blot

ELISA

Flow Cytometry

Immunoprecipitation

Protein Purification

Cell Culture

Antibody Production

Genotyping

Sequencing

Mass Spectrometry

Protein Structure

Mink Translation

Product Catalog

 Enzo Diagnostics



Contents

<i>Ordering Information</i>	<i>2</i>
<i>Company Background</i>	<i>3</i>
<i>BioProbe® Labeling Systems for Nucleic Acids</i>	<i>5</i>
<i>Hybridization and Detection Assay Systems</i>	<i>11</i>
<i>In situ Hybridization and Detection Assays</i>	<i>11</i>
<i>Identification of Human Papillomavirus</i>	<i>15</i>
<i>BioProbe® Labeled Probes</i>	<i>19</i>
<i>Dot Blot Hybridization and Detection</i>	<i>20</i>
<i>The Enzo Microplate Hybridization Assay</i>	<i>21</i>
<i>DETEK® Signal Generating Systems</i>	<i>23</i>
<i>Immunopathology Reagents</i>	<i>27</i>
<i>Monoclonal Antibodies and Detection Systems</i>	<i>27</i>
<i>Glass Fiber Filters</i>	<i>33</i>
<i>References – DNA Probe Technology</i>	<i>35</i>
<i>References – Immunopathology</i>	<i>39</i>
<i>Index</i>	<i>43</i>


The Enzo Microplate Hybridization Assays*

- ***For direct detection of HIV***
- ***For specific identification of MTB***

The Enzo Microplate Hybridization Assay System is an easy-to-use, rapid and nonradioactive kit for detecting DNA in a microtiter well. The straightforward assay procedure is based on a two probe hybridization method and is carried out entirely in the 96-well microtiter plate or, alternatively, microwell strips. The generation of a color, which can be easily measured, indicates a positive reaction.

No specialized equipment, other than a standard microplate reader, which is present in most medical research labs, is required. In fact, for plus/minus determinations the results can be read by eye. The assay format uses either a microtiter plate, where up to 96 assays (including controls) can be run, or well strips, where 8 assays per strip (including controls) can be run. The assay requires less than three hours.

Schematic Representation of the Enzo Microplate Hybridization Assay*

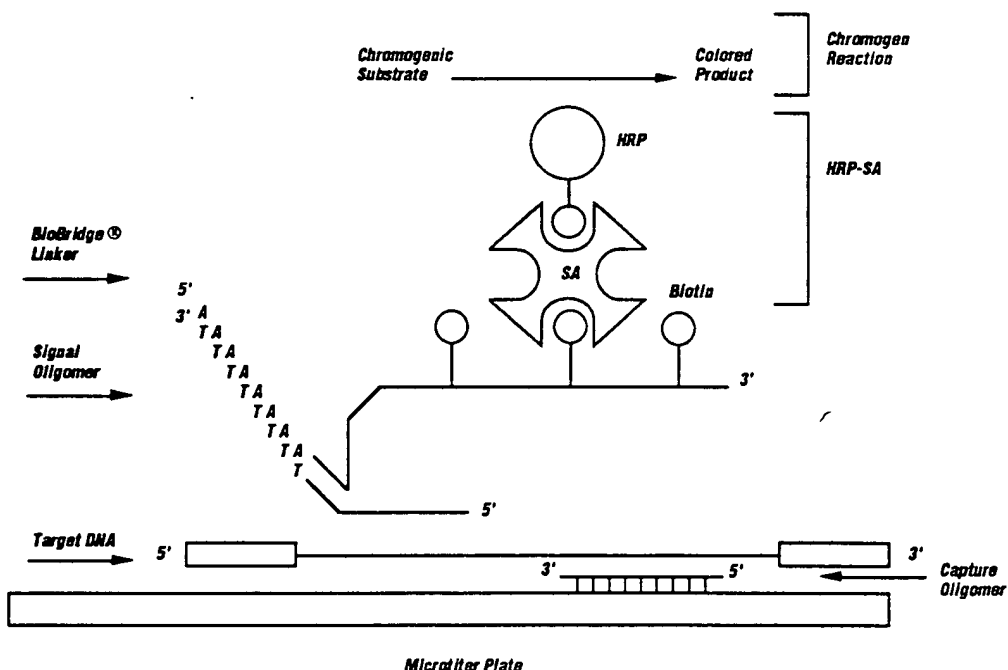

Procedure for the Enzo Microplate Hybridization Assay

This nonradioactive procedure involves hybridization of target nucleic acid to a well-bound capture probe followed by binding of a biotinylated signaling probe to the captured target. Immobilized target DNA is then visualized by reaction with a biotin-binding signal generating complex of streptavidin (SA) and horseradish peroxidase (HRP). A positive reaction is indicated by generation of color which can be measured by a microplate reader commonly used in the laboratory. Using this format, 5×10^7 to 1×10^8 copies of target sequences can be detected.

The assay utilizes several unique features. The use of probe pairs increases the specificity of the assay. Two independent hybridization events are required to generate a signal. This type

* U.S. Patent Nos. 4,711,955; 4,994,373 and Patents Pending

of assay has been found to be insensitive to the presence of cellular components other than the target DNA sequence. Since the assay is nonradioactive, the components are safe and stable, and the assay can be carried out manually or using devices developed for the automated processing of microtiter plate-based assays.

Denature Sample

Incubate sample with denaturation buffer for 15 minutes to denature the target nucleic acid sequences.

Hybridize to Well-Bound Capture Oligomer

After rinsing the wells, add sample to each well containing hybridization buffer and incubate for 90-120 minutes.

Reaction with Signal Oligomer

After the sample DNA is hybridized to the capture probe, signal probe is added and incubated for 15 minutes.

Addition of Linker

After rinsing the wells, add the linker (which supplies the biotin) to each well and incubate for 10 minutes.

Detection

Bind the streptavidin-horseradish peroxidase complex, wash the wells and then add the chromogen/substrate to generate color.

- A positive result appears as a blue color which turns to yellow upon addition of the stop solution.
- Results may be quantified by reading OD at 450.

References

1. Cook, A.F., Vuocolo, E. and Brakel, C.L. (1988) Nuc. Acid Res. 16:4077.
2. Picken, R.N., Plotch, S.J., Wang, Z., Lin, B.C., Donegan, J.J. and Yang, H.L. (1988) Molecular and Cellular Probes 2:111.

Enzo HIV Microplate Hybridization Assay*

Cat. No. 46330/plate **Using a 96 well microplate**
Cat. No. 46330/strips **Using 12 8-well strips**

The Enzo HIV Microplate Hybridization Assay Kit provides materials for the colorimetric detection of nucleic acid containing the *gag* region. HIV DNA can be assayed directly if there is sufficient target DNA present, or it can be assayed in procedures employing target amplification.

The Microplate Assay detects HIV proviral DNA with excellent sensitivity; fewer than 10 proviral sequences can be detected. Furthermore, when the assay is performed using a standard curve of HIV DNA copies, quantitative virus measurements can be done. Thus, the Microplate Assay is readily applicable to studies measuring the effect of drug treatments on virus concentration, virus concentration during the course of

infection, virus concentration in animal model studies and numerous other studies where virus quantitation is a critical parameter.

Direct and quantitative HIV detection offers a distinct advantage over serological assays. The presence of provirus DNA can be determined independent of the presence or absence of antibodies which appear several months after initial infection. Thus, this technique offers researchers the potential for identifying proviral DNA in instances where HIV positivity cannot be identified adequately by any of the current means available as, for example, in samples from HIV positive individuals who have not seroconverted.

Enzo MTB Microplate Hybridization Assay*

Cat. No. 46340/plate **Using a 96 well microplate**
Cat. No. 46340/strips **Using 12 8-well strips**

The Enzo MTB Hybridization Assay Kit provides materials for the identification of *Mycobacterium tuberculosis* complex. The assay detects only the members of the MTB complex, *i.e.*, *M. africanum*, *M. bovis*, *M. tuberculosis*, *M. microti*, but not 25 other species of *Mycobacteria*, including *M. avium*, *M. intracellulare*, *M. kansasii*, *M. chelonae*. Furthermore, 63 other bacterial species were found to be negative, including a variety of respiratory pathogens. MTB complex DNA can be assayed directly if there is sufficient target DNA present, or it can be assayed in procedures employing target amplification or culture.

While the presence of acid fast bacteria in a specimen can be highly suggestive of a mycobacterial infection, it is necessary to identify the species present. The identification of *Mycobacteria* is often more difficult than many other pathogens since the clinically relevant species do not grow rapidly in laboratory culture. Using the Microplate Assay, the researcher can specifically identify members of the MTB complex under conditions which are unrelated to the state of growth of the culture or specimen.

Oligonucleotides

Cat. No. 46331 **SK 38/SK 39**
Cat. No. 46341 **MTB 10/MTB 11**

The SK 38/SK 39 oligonucleotide pair is complementary to the HIV *gag* region. The MTB 10/MTB 11 oligonucleotide pair is complementary to sequences specific for the MTB complex.

The oligonucleotide pairs are supplied in a quantity of 2 nanomoles each oligomer.

* U.S. Patent Nos. 4,711,955; 4,994,373 and Patents Pending



Index

Cat. No.	In Situ Hybridization Assay Systems	Page
32881	BioPap™ Human Papillomavirus DNA Assay for Cervical Smears	18
32889	BioPap™ 16/18 DNA Assay for HPV Types 16/18 in Cervical Smears	18
32892	BioPap™ Human Papillomavirus DNA Typing Assay for Cervical Specimens (Types 6/11, 16/18, 31/33/51)	18
32883	BioPap™ Human Papillomavirus DNA Typing Assay Cervical Specimen Transport Kit	18
32879	PathoGene® <i>in situ</i> Screening Assay for Human Papillomavirus	13,17
32895	PathoGene® <i>in situ</i> Typing Assay for Human Papillomavirus (Types 6/11, 16/18, 31/33/51)	13,17
32877	PathoGene® DNA Probe Assay for Identification of Human Papillomavirus (Types 6/11, 16/18, 31/33/51)	13,17
32878	PathoGene® HPV 18 DNA Probe Reagent	17
32871	PathoGene® DNA Probe Assay for Identification of Adenovirus	12
32872	PathoGene® DNA Probe Assay for Identification of Cytomegalovirus	12
32873	PathoGene® DNA Probe Assay for Identification of Epstein-Barr Virus	13
32874	PathoGene® DNA Probe Assay for Identification of Hepatitis B Virus	13
32875	PathoGene® DNA Probe Assay for Identification of Herpes Simplex Virus	13
32876	PathoGene® DNA Probe Assay for Identification of <i>Chlamydia trachomatis</i>	13
32800	PathoGene® Tissue Preparation Kit	13

Cat. No.	Hybridization and Detection Assays	Page
46305	Dot Blot Hybridization and Detection Assay Kit	20
46305/C	Dot Blot Hybridization and Detection Assay Kit, Control DNA Pack	20
46307	Dot Blot Hybridization and Detection Assay Kit, CMV Control DNA Pack	20
46308	Dot Blot Hybridization and Detection Assay Kit, HBV Control DNA Pack	20
44300	Dot Blot Manifold	20
46330	HIV Microplate Hybridization Assay	22
46331	SK 38/SK 39 Oligonucleotide pair for HIV	22
46340	MTB Microplate Hybridization Assay	22
46341	MTB 10/MTB 11 Oligonucleotide pair for MTB	22

Cat. No.	Hybridization Accessories	Page
31871	Adenovirus Control Slide	14
31872	Cytomegalovirus Control Slide	14
31873	Epstein-Barr Virus Control Slide	14
31875	Herpes Simplex Virus Control Slide	14
31876	<i>Chlamydia trachomatis</i> Control Slide	14
31877	HPV 16 Probe Control Slide	14
31802	Pretreated slides, teflon coated, single well	14
31500	Heating Block for use with DNA Probe Assays 110V, 50/60HZ	13
31508	Heating Block for use with DNA Probe Assays 220V, 50HZ	13

Cat. No.	BioProbe® Labeled Probes	Page
40834	Adenovirus	19
40835	Cytomegalovirus	19
40836	Epstein-Barr Virus	19
40837	Hepatitis B Virus	19
40838	Herpes Simplex Virus	19
40839	<i>Chlamydia trachomatis</i>	19
40840	Lambda	19

Keratin-903

Melanoma

DETEK[®]

ImmunODETEK[®]

Dia 1-1 d11TD



Eizo Diagnostics, Inc.

61 Executive Drive, Farmingdale, NY 11735

(800) 224-7705, in US and Canada

(516) 496-8080, in NY State